


HACKER INSTRUMENTS

# HCM 6000 AUTOMATED GLASS COVERSLIPPER

**The HCM 6000...  
designed for carefree,  
user friendly, automated  
glass coverslipping.**

**Unique Non-Pressurized  
Media Dispensing  
Small Footprint for Space  
Savings  
Microwaveable Slide  
Baskets**

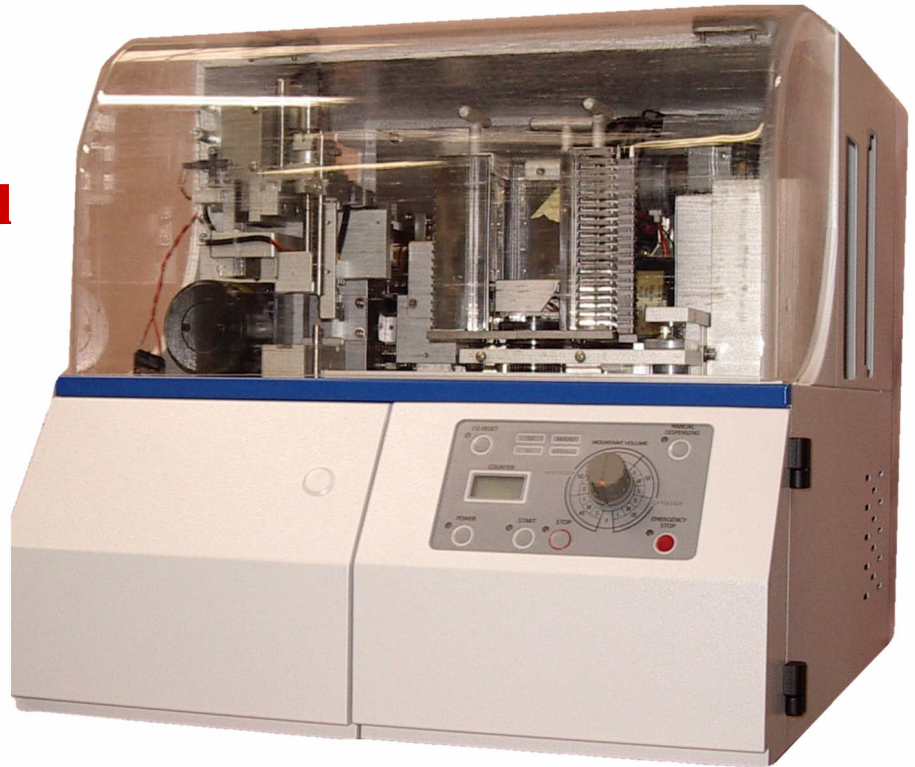

---

20 years of automated glass coverslipping experience have led to the development of the **HCM 6000 Coverslipper**. The HCM6000 offers a 60 slide loading capacity.

The unique dispensing system eliminates the formation of air bubbles under the coverglass.

All new microwaveable slide baskets interface with most automated staining machines. Available with semi-disposable, ready to load, pre-filled coverglass and mountant dispensers. Compact footprint requires minimal counterspace.

---

***Hacker Instruments and the HCM6000...  
automated, glass coverslipping you can rely on!***

## Technical Specifications - HCM6000 Glass Coverslipper

<b>SYSTEM</b>	- Control - Environment	microprocessor temperature 10 - 40° C humidity 30 - 90% RH
<b>COVERSLIPPING</b>	- Time/slide - Method	less than 14 seconds per slide Gradual Pressure
<b>GLASS SLIDE</b>	- dimensions  - Capacity - Holding Method (before coverslipping) - Storage (after coverslipping)	25 x 76 mm (1 x 3 in.) - polished edge - 0.9 - 1.2mm thickness HCM6000 up to 60 Specifically designed basket set in enclosed holding trough with xylene or xylene substitute stored in specifically designed storage rack
<b>COVERSLIP</b>	- Width - Lengths - Thickness - Capacity	24 mm (preferred) & 25mm 40mm, 50mm, 55mm, 60mm 0.13mm 200
<b>MOUNTANT</b>	- Method - Type recommended  - Viscosity - Consumption - Capacity	in-line dropping Malinol. Please check other commercial brands with manufacturer 50 to 600 cp 0.02-0.8cc adjustable 200cc internal dispensing system/500cc external dispenser
<b>FUME CONTROL</b>	- Standard - Optional	built in exhaust fan External, activated Carbon Filter Unit
<b>POWER</b>	- Standard - Optional	AC117±17V 50/60 Hz 300VA 220V 50/60 Hz 300VA
<b>ALARM</b>	- Error Message Indication for :	coverglass pick up slide pick up slide storage
<b>INDICATOR</b>	for reloading consumables & maintenance	shows quantity of slides coverslipped
<b>DIMENSIONS</b>		405 x 530 x 430 mm/ 40.5 x 53 x 43 cm
<b>WEIGHT</b>		35kg

*Call us to request a complete  
catalogue of H/I products  
for HistoPathology,  
1-800-442-2537*



**HACKER Instruments & Industries Inc.**  
1132 Kincaid Bridge Road  
Winnsboro, SC 29180 USA  
Tel: 803-712-6100 : Fax: 803-712-6116

[www.hackerinstruments.com](http://www.hackerinstruments.com)

[email: hackerlab@aol.com](mailto:hackerlab@aol.com)